

Spider IR 1 μm Central Wavelength

Towards Short Pulses at 1 μm Central Wavelength

- The Spider IR is a precision tool optimized for the complete spectral and temporal characterization of laser pulses in the infrared. Based on the patented Spider* technology, it extends the existing range of APE Spider models to cover longer pulses, between 30 and 500 fs, at a central wavelength of around 1 μm .
- It also supports detection of the chirp sign for stretched pulses greater than 2 ps width, making it a smart choice for the alignment of pulse compressors.
- With its two internal spectrometers (for fundamental spectrum and upconverted interferogram) the Spider IR is able to simultaneously measure and analyze both spectra needed for pulse reconstruction, by using the same pulse. This gives it true single-shot capability.
- Furthermore, the Spider IR control software supports real-time calculation of the temporal amplitude and phase. The user-friendly design features highly automated software to guide the operator through calibration and alignment procedures and enable measurements to be executed with a minimum of data input.



- Best choice for pulses between 30 and 500 fs at 1 μm central wavelength
- Spectral intensity and phase plus temporal intensity and phase measurement
- Real-time and true single-shot measurement of intensity and phase
- High level of automated software support and internal camera-assisted alignment
- Full software suite included

*Spectral Phase Interferometry for Direct Electric-field Reconstruction; International Patent No.: EP 1000315, WO 1999/006794

Spider IR Specifications

Specifications

Wavelength Range	970 ... 1070 nm *
Typical Application	Characterization of lasers with small bandwidths, larger pulse widths, e.g. 1 μ m laser
Spectral Bandwidth	> 6 ... 50 nm
Pulse Width	30 ... 500 fs; 5 ps for chirp direction measurement only
Laser Repetition Rate	Any; Single Shot
Input Polarization	Linear horizontal
Input Power	~ 100 mW at e.g. 80 MHz; ~ 20 mW at e.g. 1 kHz
Input Trigger	TTL for $f < 10$ Hz
Connection	USB
Software	Included; Features e.g. <ul style="list-style-type: none"> ▪ Alternative interferogram demodulation analysis: Fourier / Wavelet ▪ E-field plot ▪ Peak power calculation ▪ Measurement of phase differences (Dispersion measurement) ▪ Spectral phase derivation up to fourth order ▪ Simulation of additional theoretical dispersion (GDD, TOD, FOD) ▪ Spectrogram (X-FROG, SHG-FROG) and Wigner trace representation of the pulse

Options

- External Beam Splitter and Beam Routing Kit For lasers with low repetition rates (kHz or less)
- Pre-installed Notebook Available
- Wavelength Ranges Customized wavelength ranges on request

Dimensions

561 x 289 x 320 mm (W x H x D) (See appendix for details)

* Other wavelength ranges on request


Appendix Technical Drawings

All Dimensions in mm

Compact LX Spider

- Compact version, ideal for the characterization of Ti:Sa lasers


page 3

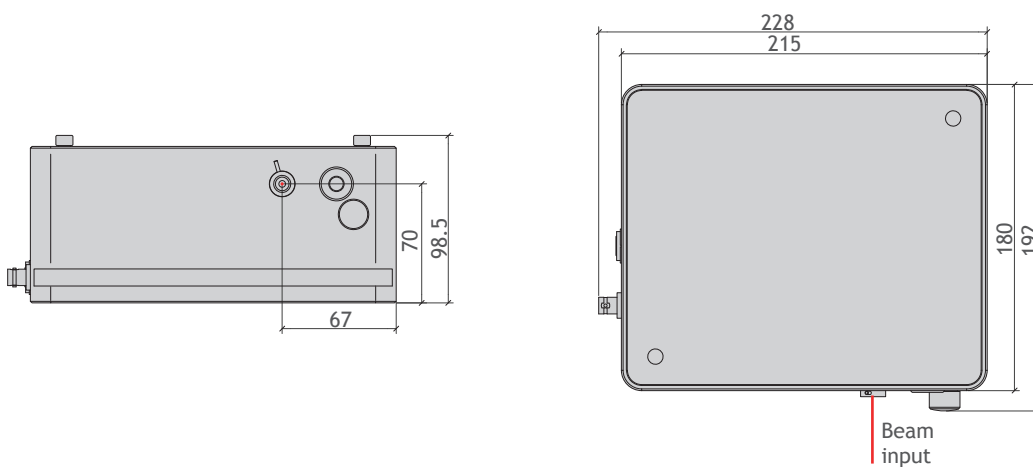

FC Spider

- Measurement of very short pulses with only a few cycles

page 5

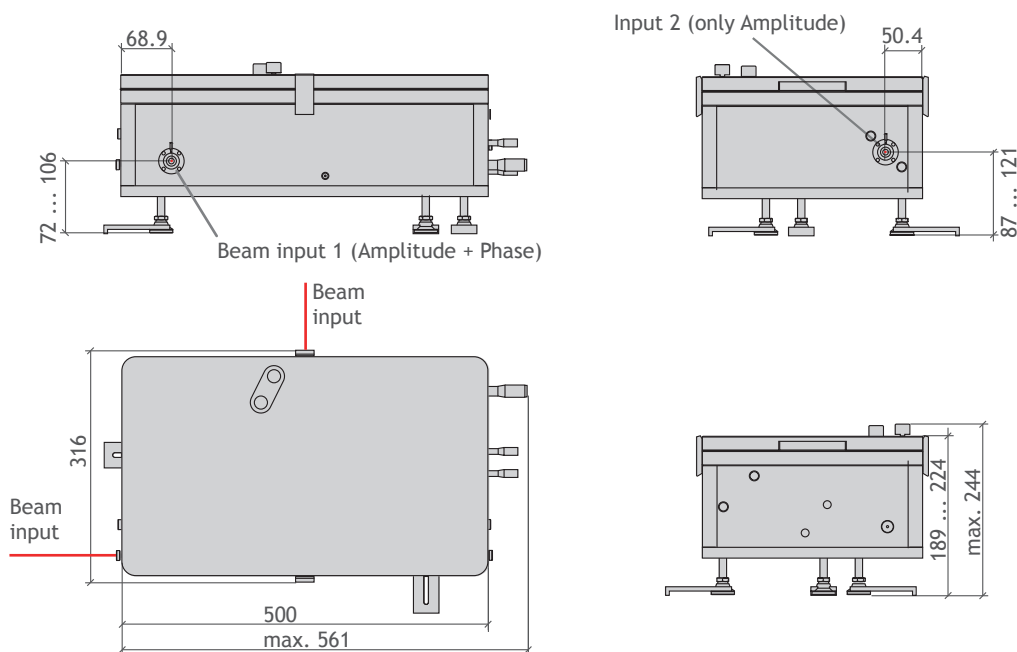


IR Spider Technical Drawings

IR Spider

- Measurement at central wavelengths of around $1 \mu\text{m}$

[page 7](#)

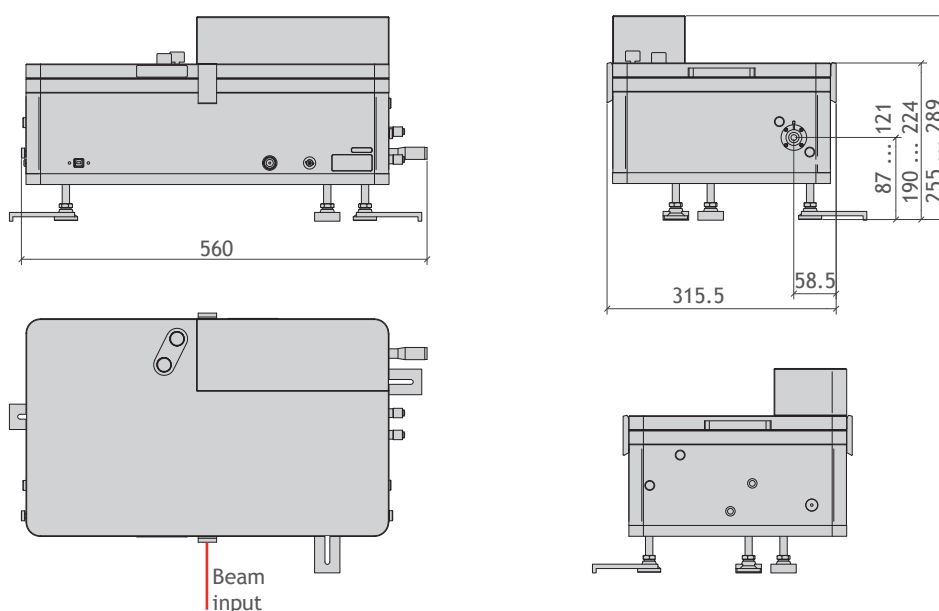

Similar Products

- pulseCheck - Autocorrelator multitalent for any task
- Mini PD - Autocorrelator routine tasks with a fixed wavelength range
- Carpe - Autocorrelator first choice for multiphoton microscopy
- waveScan - High resolution spectrometer
- peakDetect - Pulse quality monitoring

Contact

APE Angewandte Physik & Elektronik GmbH

Plauener Str. 163-165 | Haus N | 13053 Berlin | Germany
 T: +49 30 986 011-30
 F: +49 30 986 011-333
 E: sales@ape-berlin.de
 www.ape-berlin.de

PHOTO www.phototechnica.co.jp
TECHNICA フォトテクニカ株式会社
 〒336-0017 埼玉県さいたま市南区南浦和 1-2-17
 TEL:048-871-0067 FAX:048-871-0068
 e-mail:voc@phototechnica.co.jp

APE follows a policy of continued product improvement. Therefore, specifications are subject to change without notice.

© APE GmbH | November 2016 | Rev. 3.1.1