Diode Pumped Sub-Nanosecond Passively Q-Switched Laser MPL1310

FEATURES

- > Up to **1 mJ** pulse energy at **1064 nm**
- > Short pulse duration < 350 ps
- > Variable up to **100 Hz** repetition rate
- > Ultra-compact
- Passively Q-switched
- > Average power 100 mW
- > High peak power 2.8 MW
- > Guaranteed > 3 Gshot lifetime
- Other wavelengths (e.g. 1053 nm, 1342 nm, 671 nm, 447 nm) are available

MPL1310 series DPSS passively Q-switched sub-nanosecond lasers deliver high peak powers at 100 Hz repetition rate. Short laser cavity is fixed on thermo-stabilized and controlled baseplate which gives extremely stable output parameters performance. Small footprint is welcome point for integration into OEM lasers. Sub-nanosecond pulse duration of < 350 ps, high pulse energy more than 1 mJ and variable repetition rate from 1 Hz to 100 Hz covers many applications like pollution monitoring, DNA analysis, supercontinuum generation and many others. Due to short pulse duration and high pulse energy laser delivers high peak power which is up to 2.8 MW. Optional conversion to green (532 nm) and ultraviolet (355 nm, 266 nm) is also available.

APPLICATIONS

- Laser-Induced Breakdown Spectroscopy (LIBS)
- > Time resolved fluorescence measurements
- > DNA analysis
- > Pollution monitoring
- > Remote sensing
- > Supercontinuum generation
- > Ignition of gas mixtures




MODEL	MPL1310	MPL1310-M0
Pulse energy		
at 1064 nm	1 mJ	0.2 mJ
at 532 nm	0.5 mJ	0.1 mJ
at 355 nm	0.25 mJ	0.05 mJ
at 266 nm	0.15 mJ	_
Typical pulse duration	< 350 ps ²⁾	< 250 ps ²⁾
Pulse to pulse energy stability (RMS)		
at 1064 nm	< 1 % ³⁾	
at 532 nm	< 2.5 % ³⁾	
at 355 nm	< 3.5 % ³⁾	
at 266 nm	< 5.0 % ³⁾	
Power drift	± 3.0 % ⁴⁾	
Pulse repetition rate ⁵⁾	1 – 100 Hz	
Beam profile	Close to Gaussian	
Beam divergence 6)	< 6 mrad	
Polarization	Linear, horizontal at 1064 nm	
Spectral linewidth	SLM	
Beam pointing stability 7)	< 40 µrad	
Typical beam diameter ⁸⁾	1.5 mm	1 mm
Jitter	~ 2 µs RMS ⁹⁾	


DIMENSIONS

Laser head (W×L×H)	113 × 162.5 × 45.5 mm
Controller unit (W×L×H)	75 × 200 × 70 mm (OEM version)
Power adapter (W×L×H)	119 × 179 × 30 mm (OEM version)

OPERATING REQUIREMENTS

Cooling requirements	Air cooled	
Ambient temperature	15 – 30 °C	
Relative humidity (non-condensing)	10 - 80 %	
Mains voltage	90 – 230 VAC, single phase, 47 – 63 Hz ¹⁰⁾	
Power consumption	< 20 W	

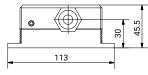

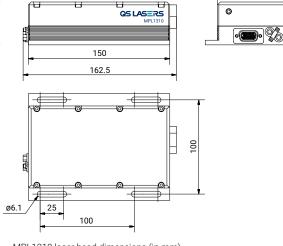

www.phototechnica.co.j

りォトテクニカ株式会社

87


MPL1310 laser head dimensions (in mm)

+370 620 42058 @ sales@qslasers.com
 Savanoriu Ave. 231, Vilnius LT-02300, Lithuania



- ¹⁾ Due to continuous improvements all specifications are subject to change. Unless stated otherwise all specifications are measured at 1064 nm.
- 2) FWHM level at 1064 nm. Other pulse duration is available by request. Please inquire for detailed specifications.
- ³⁾ Averaged from 60 seconds time interval in 5 series.
- ⁴⁾ Over 8-hour period after max 5 minutes of warm-up when ambient temperature variation is less than ±2 °C.
- ⁵⁾ Factory-set pulse repetition rate is fixed at 100 Hz repetition rate. Higher repetition rates are available, please inquire for more details.
- ⁶⁾ Full angle measured at the 1/e² level.
- 7) RMS value measured from 1000 shots.
- ⁸⁾ Beam diameter is measured 20 cm from laser output at the 1/e² level.
- ⁹⁾ In respect to Q-switch triggering rising edge pulse.
- Laser can be powered from appropriate 12 VDC power source. Inquire for details.



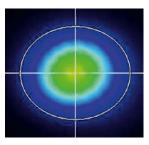

Typical beam intensity profile (20 cm from laser output) of MPL1310 series lasers

Diode Pumped Sub-Nanosecond Passively Q-Switched Laser MPL1510

FEATURES

- > Up to **1 mJ** pulse energy at **1064 nm**
- > Short pulse duration < 500 ps
- > Variable up to **100 Hz** repetition rate
- > Ultra-compact
- Passively Q-Switched
- > Average power 100 mW
- > High peak power **2 MW**
- > Guaranteed > 3 Gshot lifetime
- Other wavelengths (e.g. 1053 nm, 1342 nm, 671 nm, 447 nm) are available

APPLICATIONS

- Laser induced breakdown spectroscopy (LIBS)
- > Time resolved fluorescence measurements
- > DNA analysis
- > Pollution monitoring
- > Remote sensing
- > Supercontinuum generation
- > Ignition of gas mixtures

MPL1510 series DPSS passively Q-switched sub-nanosecond lasers deliver high peak powers at 100 Hz repetition rate. Short laser cavity is fixed on thermo-stabilized and controlled baseplate which gives extremely stable output parameters performance. Small footprint is welcome point for integration into OEM lasers. Sub-nanosecond pulse duration of < 500 ps, high pulse energy more than 1 mJ and variable repetition rate from 1 Hz to 100 Hz covers many applications like pollution monitoring, DNA analysis, supercontinuum generation and many others. Due to short pulse duration and high pulse energy laser delivers high peak power which is up to 2.8 MW. Optional conversion to green (532 nm) and ultraviolet (355 nm, 266 nm) is also available.






MODEL	MPL1510	
Pulse energy		
at 1064 nm	1 mJ	
at 532 nm	0.5 mJ	
at 355 nm	0.25 mJ	
at 266 nm	0.15 mJ	
Typical pulse duration	< 500 ps ²⁾	
Pulse to pulse energy stability (RMS)		
at 1064 nm	< 1 % ³⁾	
at 532 nm	< 2.5 % ³⁾	
at 355 nm	< 3.5 % ³⁾	
at 266 nm	< 5.0 % ³⁾	
Power drift	± 3.0 % ⁴⁾	
Pulse repetition rate ⁵⁾	1 – 100 Hz	
Beam profile	Close to Gaussian	
Beam divergence 6)	< 6 mrad	
Polarization	Linear, horizontal at 1064 nm	
Spectral linewidth	SLM	
Beam pointing stability 7)	< 40 µrad	
Typical beam diameter ⁸⁾	2 mm	
Jitter	~ 2 µs RMS ⁹⁾	


DIMENSIONS

Laser head (W×L×H)	113 × 162.5 × 45.5 mm
Controller unit (W×L×H)	75 × 200 × 70 mm (0EM version)
Power adapter (W×L×H)	119 × 179 × 30 mm (OEM version)

OPERATING REQUIREMENTS

Air cooled
15 – 30 °C
10 - 80 %
90 – 230 VAC, single phase, 47-63 Hz ¹⁰⁾
< 20 W

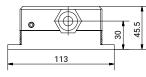

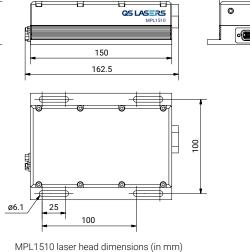

87

D• %

www.phototechnica.co.jp

フォトテクニカ株式会社




😂 +370 620 42058 @ sales@qslasers.com 🌐 www.qslasers.com

💡 Savanoriu Ave. 231, Vilnius LT-02300, Lithuania

¹⁾ Due to continuous improvements all specifications are subject to change. Unless stated otherwise all specifications are measured at 1064 nm.

- PWHM level at 1064 nm. Shorter pulse duration (< 350 ps) is available by request. Please inquire for detailed specifications.
- Averaged from 60 seconds time interval in 5 series.
- ⁴⁾ Over 8-hour period after max 5 minutes of warm-up when ambient temperature variation is less than ±2 °C.
- ⁵⁾ Factory-set pulse repetition rate is fixed at 100 Hz repetition rate. Higher repetition rates are available, please inquire for more details.
- ⁶⁾ Full angle measured at the 1/e² level.
- ⁷⁾ RMS value measured from 1000 shots.
- ⁸⁾ Beam diameter is measured 20 cm from laser output at the 1/e² level.
- ⁹⁾ In respect to Q-switch triggering rising edge pulse.
- Laser can be powered from appropriate 12 VDC power source. Inquire for details.



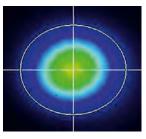

Typical beam intensity profile (20 cm from laser output) of MPL1510 series lasers



Diode Pumped Picosecond Passively Q-Switched Lasers MPL2310

FEATURES

- > More than 2 mJ pulse energy at 1064 nm
- Short pulse duration < 500 ps</p>
- > 1 100 Hz repetition rate
- Ultra-compact
- Passively Q-switched
- > Average power 200 mW
- > High peak power > 4 MW
- > Guaranteed > 3 Gshot lifetime
- Other wavelengths (e.g. 532 nm, 355 nm, 266 nm) are available

MPL2310 series DPSS passively Q-switched sub-picosecond lasers deliver high peak powers > 5 MW at 100 Hz repetition rate. Short laser cavity is fixed on thermostabilized and controlled baseplate which gives extremely stable output parameters performance. Small footprint is welcome point for integration into OEM lasers. Subnanosecond pulse duration of < 350 ps, high pulse energy more than 2 mJ, variable repetition rate from 1 Hz to 100 Hz covers many applications like pollution monitoring, DNA analysis, supercontinuum generation and many others.

Due to short pulse duration and high pulse energy laser delivers high peak power which is up to 5 MW. Optional conversion to green (532 nm) and ultraviolet (355 nm, 266 nm) is also available.

APPLICATIONS

- Laser-induced breakdown spectroscopy (LIBS)
- > Time resolved fluorescence measurements
- > DNA analysis
- > Pollution monitoring
- > Remote sensing
- > Supercontinuum generation
- > Ignition of gas mixtures







MODEL		MDI 1010 / MDI 1510	
MODEL	MPL2310 / MPL2510	MPL1310 / MPL1510	
Pulse energy			
at 1064 nm	2 mJ	1 mJ	
at 532 nm	1 mJ	0.5 mJ	
at 355 nm	0.5 mJ	0.25 mJ	
at 266 nm	0.25 mJ	0.15 mJ	
Typical pulse duration	<350 ps/ <500 ps ²		
Pulse to pulse energy stabilit	y (RMS)		
at 1064 nm	< 1.0 % ³⁾		
at 532 nm	< 2.0 % ³⁾		
at 355 nm	< 3.0 % ³⁾		
at 266 nm	< 4.0 % ³⁾		
Power drift	± 3.0 % ⁴⁾		
Pulse repetition rate ⁵⁾	1 – 100 Hz		
Beam profile	close to Gaussian		
Beam divergence 6)	< 6 mrad		
Polarization	linear, horizontal at 1064 nm		
Spectral linewidth	SLM		
Beam pointing stability 7)	< 10 µrad		
Typical beam diameter ⁸⁾	1.5 mm		
Optical jitter	~ 2 µs RMS ⁹⁾		
DIMENSIONS			
	125 × 295 × 76 mm (with harmonics)		
Laser head (W×L×H)	99 × 174 × 45.5 mm (OEM version)		
	257 × 271 × 153 mm		
Controller unit (W×L×H)	75 × 200 × 70 mm (OEM version)		


OPERATING REQUIREMENTS

Cooling requirements	air cooled	
Ambient temperature	15 – 30 °C	
Relative humidity	10 - 80 % (non-condensing)	
Mains voltage	100 – 230 VAC, single phase, 50 – 60 Hz ¹⁰⁾	
Power consumption	< 20 W	< 10 W


- ¹⁾ Due to continuous improvements all specifications are subject to change. Unless stated otherwise all specifications are measured at 1064 nm.
- ²⁾ FWHM level at 1064 nm.
 ³⁾ Averaged from 60 seconds time interval in 5 series.
- ⁴⁾ Over 8-hour period after max 5 minutes of warm-up when ambient temperature variation is less than ±2 °C.
- Factory-set pulse repetition rate is fixed at 100 Hz repetition rate. Higher repetition rates are available, please inquire for more details.
- ⁶⁾ Full angle measured at the 1/e² level. Lower beam divergence is available upon request, please inquire for more details.
- 7) RMS value measured from 1000 shots.
- ⁸⁾ Beam diameter is measured 20 cm from laser output at the 1/e² level.
- ⁹⁾ In respect to Q-switch triggering rising edge pulse.

¹⁰⁾ Laser can be powered from appropriate 12 VDC power source. Inquire for details.




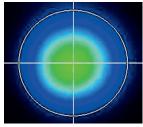

Typical beam intensity profile (20 cm from laser output) of MPL2310 series lasers

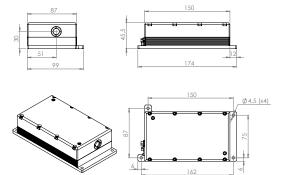

MPL2310 series laser head dimensions OEM version (in mm)

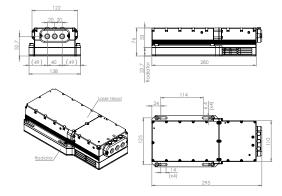

MPL2310 series laser head dimensions (in mm)

www.phototechnica.co.jp

フォトテクニカ株式会社

+370 620 65424 @ sales@qslasers.com
 Waujoji Riovoniu str. 25B, Vilnius LT-03153, Lithuania



Diode Pumped Sub-Nanosecond Actively Q-Switched Laser MPL15100

FEATURES

- > More than 0.5 mJ pulse energy at 1064 nm
- Short pulse duration < 700 ps</p>
- > 1000 Hz repetition rate
- 532 nm, 355nm, 266 nm wavelengths as standard option
- > Actively Q-Switched
- > High Peak Power 0.7 MW
- Other wavelengths (e.g. 1053 nm, 1342 nm, 671 nm, 447 nm) are available

MPL15100 series robust DPSS actively O-switched sub-nanosecond lasers deliver multi-kW peak powers, less than 1 ns pulse duration at 1000 Hz repetition rate. Short innovative laser cavity with is fixed on thermostabilized baseplate which gives extremely stable output parameters performance. Small footprint is welcome point for integration into OEM lasers. Sub-nanosecond pulse duration of < 700 ps with near transform limited spectral linewidth at repetition rates up to 1 kHz with low timing jitter of <200 ps and energies more than 500 µJ covers broad spectrum of applications starting from LIBS, laser induced fluorescence to many others. Standard optional harmonics generator to green (532 nm) and ultraviolet (355 nm, 266 nm) is also available.

APPLICATIONS

- > OLED repair
- Marking
- > Nonlinear Optics
- > Seeding laser amplifiers
- > Pollution Monitoring
- Remote sensing





MODEL	MPL15100	MPL15100-1K	
Pulse energy:	'	1	
at 1064 nm	0.5	mJ	
at 532 nm	0.25	5 mJ	
at 355 nm	0.15	5 mJ	
at 266 nm	0.05	5 mJ	
Typical pulse duration	< 700 ps		
Pulse to pulse energy stabilit	y (RMS):		
at 1064 nm	< 0.5 % ²⁾		
at 532 nm	< 1.0 % ²⁾		
at 355 nm	< 1.5 % ²⁾		
at 266 nm	< 2.0 % ²⁾		
Typical pulse duration	≤ 700 ps ³⁾		
Power drift	± 3.0 % ⁴⁾		
Pulse repetition rate ⁵⁾	100 Hz	1000 Hz	
Beam spatial profile	Close to Gaussian		
Beam divergence 6)	< 4 mrad		
Polarization	Linear, horizontal at 1064 nm		
Spectral linewidth	SLM		
Beam pointing stability 7)	< 50 µrad		
Typical beam diameter ⁸⁾	1.2 mm		
Optical jitter	< 0.3 ns ⁹⁾		
DIMENSIONS			


- ¹⁾ Due to continuous improvements all specifications are subject to change. Unless stated otherwise all specifications are measured at 1064 nm.
- ²⁾ Averaged from 60 seconds time interval.
- FWHM level at 1064 nm.
 Over 8-hour period after max 5 minutes of warm-up when ambient temperature variation is less than ±2 °C.
- Factory-set pulse repetition rate is fixed at max repetition rate. Higher repetition rates are available, please inquire for details.
- ⁶⁾ Full angle measured at the 1/e² level.
- 7) RMS value measured from 1000 shots.
- 8) Beam diameter is measured 20 cm from laser output at the 1/e² level.
- ⁹⁾ In respect to Q-switch triggering rising edge pulse.




Laser head (W×L×H)	163 × 295 × 53.5 mm
Controller unit (W×L×H)	257 × 271 × 153 mm
Cable cord length	1 m

OPERATING REQUIREMENTS

Cooling requirements	air cooled	
Ambient temperature	15 – 30 °C	
Relative humidity	10 - 80 % (non-condensing)	
Mains voltage	100 – 240 VAC, single phase, 50 – 60 Hz	
Power consumption	< 10 W	< 100 W

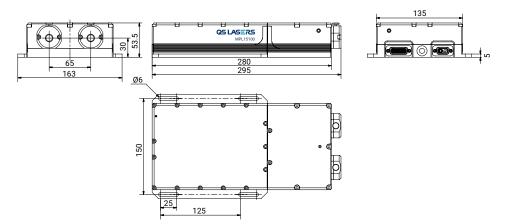


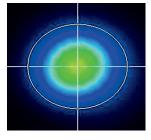

Typical beam intensity profile (20 cm from laser output) of MPL15100 series lasers

MPL15100 laser head dimensions with attached harmonics unit (in mm)

+370 620 42058 @ sales@qslasers.com
 Savanoriu Ave. 231, Vilnius LT-02300, Lithuania

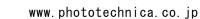

フォトテクニカ株式会社

